



Jorgensen Laboratories

Happy New Year!!

from

Jorgensen Laboratories

January 2009

Volume 14

Issue 1

[Page 1](#)

Specials

- Annual Suture Sale

2009 Catalog

[Page 2](#)

Upcoming Special

- Cal-Gel, Kato-Gel

- Pet Pillers

[Page 3](#)

Clinic Highlight

New Products

- Infusion Screw

[Page 4](#)

New Products

- Canine Ventilator Circuit

- Absorbent Kennel Pads

- BioCam Video Otoscope

[Page 5](#)

New Products

- Enteral Feeding Tubes

- Thora Cath

[Page 6](#)

New Products

- Round Bottom Tubes

- Equine Safety Restraint

[Page 7](#)

New Products

- Large Bore IV Set

[Page 8](#)

- Company Information

JorVet's Annual Suture Sale


The new year just would not be ushered in right without the JorVet annual J-91 Braun suture sale. The sale runs January 1-31, 2009. JorVet will supply one cassette free with every ten ordered (must be the same kind and size). This includes the popular Braunamid, Catgut, and Stainless Steel suture.



Needle suture will be 10% off the regular prices.

JorVet Catalog for 2009

Available in color

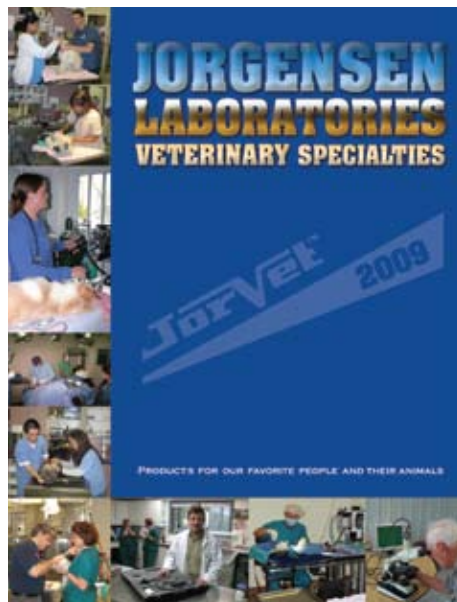

The start of a new year always means the release of a new catalog. Packed with over 360 pages of JorVet merchandise, the catalog has over 20 pages of new and innovative products for 2009.

Call for your 2009 Catalog today!!

Coming in February

Get Gross with JorVet J-294 and J-306


During the month of February, when ordering a gross of Cal-Gel (Respond) or Keto-Gel, you will receive a 10% discount. One dozen per box.



J-294

Retail

J-294 \$65.55

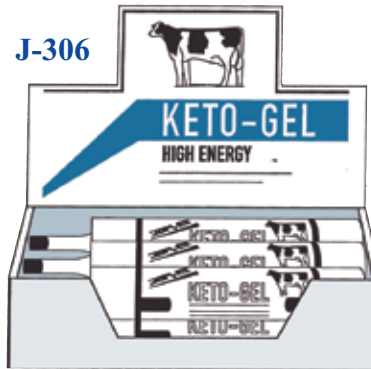

J-306

Retail

J-306 \$78.00



Pilling Made Easy with JorVet

Also during February, the balling gun style JorVet Pet Piller will be 10% off.



J-384

Retail

J-384 \$16.80
per dozen

J-384x Buster pet pillers

new and improved model

- soft rubber jaws to reduce the chance for injury and are more comfortable for the patient
- finger rings for better control
- can draw water into pet piller to aid swallowing
- attractive display boxes hold 10 pet pillers
- transparent barrel

J-384x

Retail
\$15.00



Prices listed do not reflect the discount.

Clinic Highlight

Aspen Meadow Veterinary Specialists

Dr. Matt Rooney, Owner

Mission: To be the veterinary specialty and emergency center of choice by providing the highest quality, most advanced and compassionate veterinary care and customer service.

Goal: To enhance the lives of our patients and their human families by delivering excellent veterinary care focused on quality service.

Vision: To be the preferred veterinary specialty and emergency hospital based on our four guiding principals: Compassionate Expert Care, Teamwork, Communication, and Environmental & Social Responsibility.

Recent JorVet Purchases: Doppler, Vacuum System, Liquid Oxygen/H-tank Manifold, Oxygen Drops, Orthopedic Equipment, Oval IV Track.



Rehab Room



Surgery Suite



Liquid Oxygen Manifold & H-Tank



O2/Vacuum Outlets

Aspen Meadow Veterinary Specialists

104 South Main Street, Longmont, CO 80501 www.AspenMeadowVet.com Info@AspenMeadowVet.com

New Products

Cannulated Infusion Screw

Where infection has penetrated and lodged in the medullary cavity of bone in the distal limb it may be necessary to infuse antibiotics directly to achieve therapeutic levels where they are needed.

Drill and tap 4.5 threaded hole (J-372e plus J-409e)

Insert Cannulated Infusion screw

Attach syringe and luer connector and infuse

Screw may be left in situ for repeated administration of antibiotics.

J343 Cannulated Luer Lock Infusion Screw with Insertion Key in SS case.

Retail \$145.00



New Products

Canine Ventilator Circuit

This circuit was designed with the help of a prominent respiratory disease specialist. The circuit was designed to work in conjunction with the J-949m Nebulizer pump. Its function is to increase the inhalation percentage and decrease the exhalation of vented product. This is accomplished by: two one-way valves that restrict flow direction and a bacterial filter that traps wasted nebulized drugs. The operator will have less exposure to nebulized drugs.

J-942 Retail
\$19.50



Absorbent Kennel Pads

Not all kennel pads are created equal!

The JorVet Kennel pads are highly absorbent with a special polymer that is designed to absorb 1500ml of water. Make sure you are using a kennel pad that actually keeps the patient dry and does not have to be changed with every drop of fluid. Also useful for procedures like dental prophylaxis to help maintain a dry field. Blue color with quilted pad.



size: 23"x36"
weight 85gm
50pk

J-1082 Retail
\$15.00

BioCam 2.0™ Video Otoscope

The portability and high-resolution imaging and video recording make clinical diagnosis simple and practical. It educates the client in understanding the process resulting in increased compliance. This one of a kind video otoscope has multi-use capabilities in one device. Designed with precision and cutting edge technology.

Features:

Portable -- This device can be plugged into any Microsoft platform computer with a USB2.0 port.

Power -- No batteries or external power source is needed. It is powered through the USB port.

Light Source -- 6 LED bulbs allow 10,000 working hours



J-1100 Retail
\$850


New Products

Enteral Feeding Products

This line of products is especially designed for mid-long term enteral feeding. Features polyurethane feeding tubes.

- More durable and longer lasting than PVC tubes
- Weighted tip with metal tungsten guide end for ease of placement and retention
- Stylets included on all sizes
- Radiopaque and non-sterile
- Many of these tubes will move by peristaltic action to the jejunum after placement in the stomach, which requires longer lengths than standard feeding tubes (i.e. cats need 36" (91cm)L)
- Numerical graduations on side of the tube
- Complete directions

		Retail
J-350a	5fr x 22"	\$15.00
J-350b	6fr x 22"	\$22.50
J-350c	6fr x 36"	\$22.50
J-350d	8fr x 22"	\$22.50
J-350e	8fr x 36"	\$22.50
J-350f	8fr x 43"	\$24.00
J-350g	10fr x 36"	\$24.00
J-350h	10fr x 43"	\$25.50
J-350j	10fr x 55"	\$27.00
J-350k	12fr x 43"	\$25.50

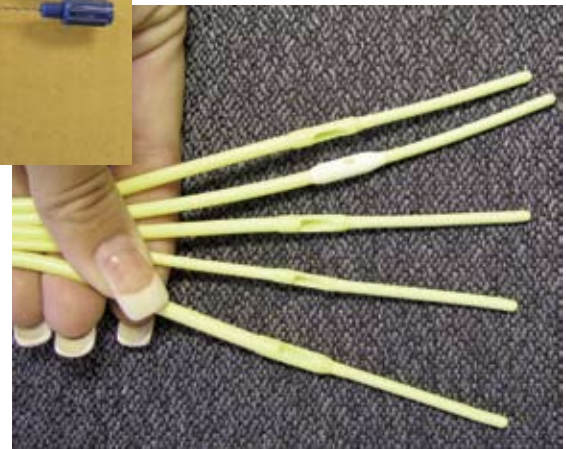


Thora-Cath® Chest Tubes

A soft silastic chest tube for pneumothorax or hemothorax. This catheter will provide drainage of the pleural space. Silicon is definitely the preferred material over pvc or nylon.

- Biocompatible silicon reduces clotting and encrustation
- Beveled end open end with two lateral eyes 1cm apart.
- Graduated at 1cm lengths
- Radiopaque stripe
- Universal hubs adapt to connectors like JorVet J282d1 Barbed reducing connector and suture wings. (Note: JorVet # J282d3 is too small.) Also suggest use of 3-way stopcock (J-462)
- Provided with stainless steel trocar
- 14cm (5½") length

		Retail
J341	9.6fr (3.2mm)	\$29.90
J341a	12fr (4mm)	\$29.90
J341b	14fr (4.67mm)	\$29.90

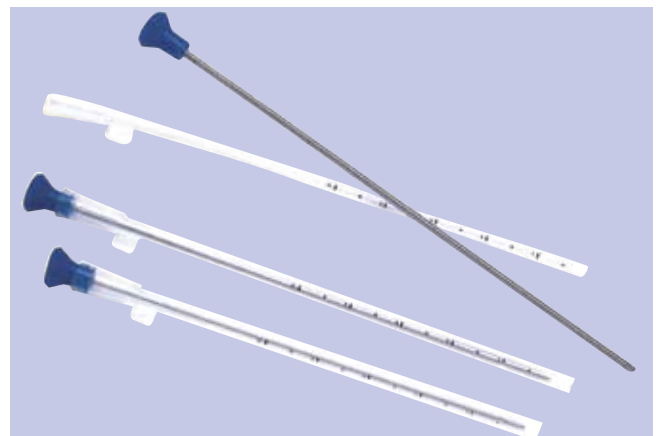

New Products

Round Bottom Centrifugal Test Tubes

- original classic design features straight edges and no cone shape end
- polystyrene -- approximately 15ml capacity
- 100mm length x 17mm diameter
- 100pk

J-546e

Retail
\$12.60



Equine Safety Restraint Straps

These straps were specifically designed to prevent/minimize any injuries to horses that are placed in equine stocks (also injuries to hospital personnel). Anyone that uses equine stocks appreciates their value but realizes that some horses may react violently to their use. They can attempt to jump over the front door or rear over backwards injuring themselves and hospital staff.

The Equine Safety Restraint Straps when placed properly over the horses shoulders will prevent lunging forward and any rearing while in the stock.

A quick release feature will immediately open to remove the straps in an emergency, even under a full load.

Materials:

- Safety strap 2" nylon webbing with strength of 2,500 lbs
- Heavy metal "O" rings
- Thick 20" wide Velcro
- Quick release

J-954



New Products

Large Bore IV Administration Set


Larger animals such as horses and cattle require tremendous amounts of IV fluids in a short period of time. The standard smaller human IV sets cannot meet the large IV fluid demands of a critically ill 1,200lb horse.

The JorVet Large Bore IV Set can increase the flow rate up to almost 20-liters per hour, a necessary volume in severe colic or equine critical care emergencies.

Poiseville's Law of Fluids: The volume of a fluid through a tube varies as the fourth power of the diameter of the tubing. $V=D^4$ other things being equal.

The diameter of the large bore IV set is two times that of a conventional IV set. This translates into 16 times the flow rate.

J-400L Large IV Set. Complete. Sterile. Retail
\$30.00

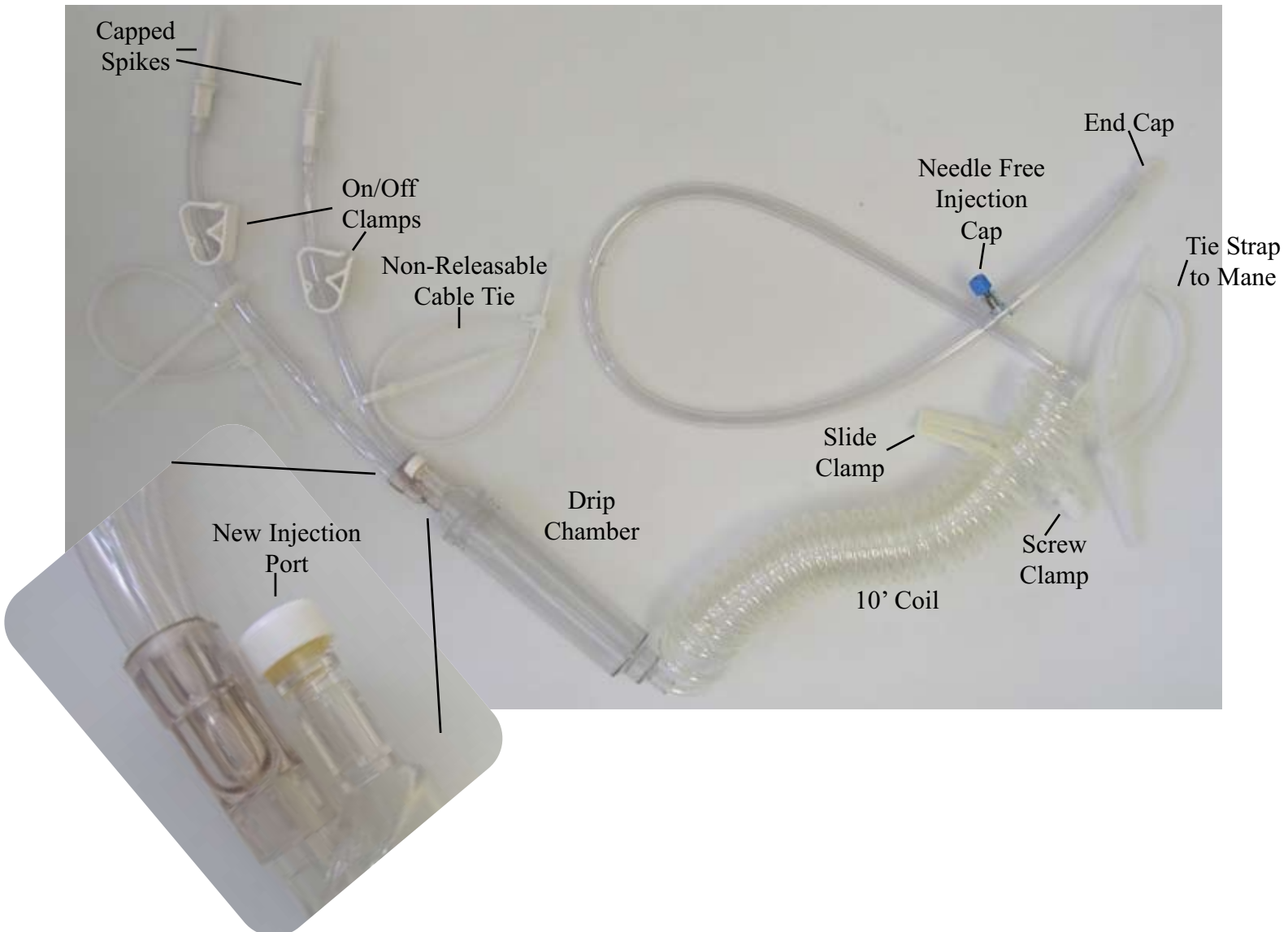



Jorgensen Laboratories, Inc.

1450 North Van Buren Avenue
Loveland, CO 80538

PRSR STD
US POSTAGE
PAID
LOVELAND, CO
80538
PERMIT NO. 26



Jorgensen Laboratories, Inc.

(800) 525-5614 * (970) 669-2500 * FAX (970) 663-5042

e-mail: Info@JorVet.com

www.JorVet.com

