


DENTAL INSTRUMENTS FOR RABBITS, GUINEA PIGS AND CHINCHILLAS



RODENT MOUTH GAG, SPATULA & CHEEK DILATOR

Without gags and dilators examination and treatment of rabbit molars is virtually impossible. Use the spatula to deflect tongue during examination of mouth. The cheek dilators allow visualisation of molars and protects cheeks from instruments. The small dilator is also useful as a small gag for guinea pigs and chinchillas.

J51rg	Rodent mouth gag	\$96.00
J51rs	Spatula	\$24.00
J51rdl	Cheek dilator large	\$21.00
J51rdx	Dilator large - long reach	\$21.00
J51rds	Cheek dilator small	\$21.00

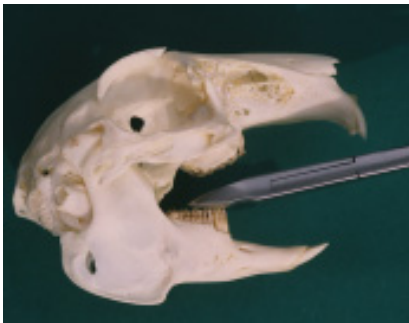

MOLAR CUTTERS

Use to debulk or for trimming spurs and hooks. The shape of the cutters is designed to allow visualisation of the teeth during cutting

J41mc	Molar cutters	\$72.00
--------------	----------------------	----------------




DIAMOND MOLAR RASP

Unsuitable for debulking but good for tidying. The diamond rasp's surface is concave to fit over the molars without slipping off. Diamond grit is the ideal material to remove tooth spurs without damaging the delicate root structure.

J41rf	Molar rasp fine - diamond	\$66.00
--------------	----------------------------------	----------------

Also available

J41gf	Molar rasp fine (guinea pigs & chinchillas).	\$40.00
J41r	Molar rasp (rabbits).	\$40.00

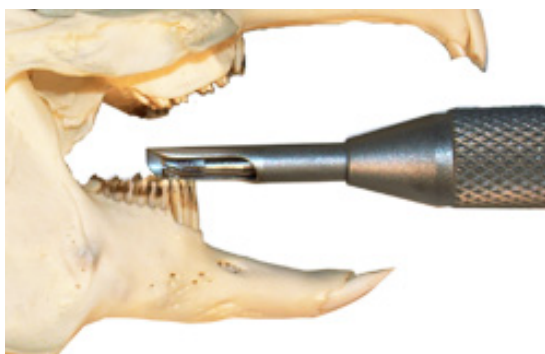

SOFT TISSUE PROTECTOR AND RABBIT BURS

J41st	Soft tissue protector (fits straight nose cone - not included).	\$106.00
J41cc	Extra long cross cut bur - for use with above.	\$8.90
J41rk	Rabbit bur kit in block.	\$65.00



DIAMOND DISC & SAFETY SHIELD

Diamond disc for trimming incisors . The safety shield is an important piece of equipment protecting both the user and the patient.

J452d7	Diamond disc	\$36.00
J41SS	Safety shield	\$99.00



CROSSLEY INCISOR LUXATOR

Designed by David Crossley to extract the very large rabbit incisors. Removal avoids the requirement for repeated cutting of overgrown mal-occluded incisors.

J41cr	Crossley incisor luxator	\$72.00
--------------	---------------------------------	----------------



CROSSLEY MOLAR LUXATOR

Designed to facilitate extraction of molars with a minimum of damage to the alveolus and surrounding soft tissue. One cutting blade is for insertion down buccal or lingual surface. The other end of the instrument is set at right angles to penetrate between molars.


J41mr	Crossley molar luxator	\$72.00
--------------	-------------------------------	----------------



MOLAR EXTRACTION FORCEPS

There is very little space inside a rabbit mouth. These forceps are very slim and angled to fit in the mouth without totally obscuring the teeth.

J41me	Molar extraction forceps	\$36.00
--------------	---------------------------------	----------------

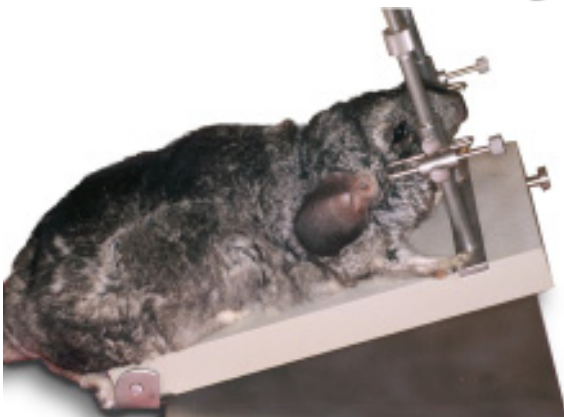

RABBIT/CHINCHILLA TABLE TOP STAND

Allows mouth to be held open in an appropriate position without assistance. The platform tips allowing best positioning for surgeon and available lighting. It works even if incisors have been removed.

J931	Table top gag.	\$395
-------------	-----------------------	--------------

ANTIBIOTIC IMPREGNATED BEADS (AIPMMA)

Make your own using Jorvet bone cement which is available with and without gentamycin. Each pack comes with instructions on how to make over 100 beads.

J913	Low viscosity cement, 40g.	\$72.00
J913g	Low viscosity cement with gentamycin, 40g.	\$72.00